



SENSORS AND AUTOMATION

Sarvesh Anand Mankar
142203013
SY Comp Div-2, S5 Batch

Practical-4: Characterization of pH, Conductivity

Aim:

1. Study the working principle of pH and conductivity sensors.
2. Calibrate the pH sensor.
3. Study the effect of temperature on pH measurement.
4. Study effect of temperature and effect of contamination on conductivity measurement.

Theory:

pH Measurement:

pH is defined as the negative logarithm of the hydrogen ion concentration. $\text{pH} = -\log(\text{H}^+)$

pH sensor:

pH is one of the most common analyses used in Process Industry. pH is actually a measurement of the activity of hydrogen ions in the sample. pH measurements run on a scale from 0 - 14, with 7.0 considered as neutral. Solutions with a pH value below 7.0 are considered as acids, and above 7.0 are designated as bases. The pH scale is logarithmic, so one unit change in pH value actually reflects a ten-fold change in the acidity.

Working principle

The pH measurement loop can be considered as a battery where the positive terminal is the measuring electrode and the negative terminal is the reference electrode.

The measuring electrode', which is sensitive to the hydrogen ions, develop a potential (voltage) directly related to the hydrogen ion concentration of the solution. The reference electrode is stable regardless of any change in the hydrogen ion concentration.

The pH meter consists of **three major components**: pH probe, Temperature probe and the meter

The pH probe consists of a glass, hydrogen-ion selective electrode, and a reference electrode, combined into a single unit. The glass electrode is specially treated for measuring the hydrogen ions, while the reference electrode is surrounded by silver chloride. It provides a "zero" or reference point for the measurement. This "zero" point means any change in potential measured at the glass electrode is attributed to hydrogen ions, and is expressed as pH. When the temperature and pH probes are immersed in the sample, the meter measures the potential difference between the glass electrode and the reference electrode. This electronic measurement is converted from millivolts to pH units, and the result appears on the display.


If the temperature probe is not used during the pH measurement, the meter will assume a temperature of 25°C.

Calibration of pH probe with buffer solution:

The calibration must be performed for buffers with pH 4.0, 7.00 and 9.2. At least a 2-point calibration must be performed at room temperature using buffers that meet the expected pH value of the sample.

If a one-point calibration is performed, measurement errors are more for the sample that is being measured. If the temperature probe is not used at the time of calibration, default value considered is 25°C.

Pretest:



Virtual
Labs
An MOE, Govt of India Initiative

HOMEPARTNERSCONTACT

Electrical Engineering > Sensors Modeling & Simulation Lab > Experiments

Aim

Theory

Pretest

Procedure

Simulation

Posttest

References

Feedback

Simulate the performance of a chemical sensor

When the hydrogen ion concentration of a solution increases the solution is called as

- ☒ a: Acidic
- ☐ b: Basic
- ☐ c: Neural
- ☐ d: Dialute

The pH of a neutral solution is

- ☐ a: 0
- ☐ b: 10
- ☐ c: 14
- ☒ d: 7

The pH scale generally ranges from

- ☒ a: 0 to 14
- ☐ b: -14 to 14
- ☐ c: -7 to 7
- ☐ d: -10 to 10

PH of distilled water is

- ☐ a: 0
- ☒ b: 7
- ☐ c: 14
- ☐ d: 10

Buffer solutions are used as a means of keeping pH at a nearly _____ value in a wide variety of chemical applications.

- ☐ a: High
- ☒ b: Constant
- ☐ c: Neutral
- ☐ d: Flexible

Conductivity of a pure and distilled water is

- ☐ a: High
- ☐ b: Low
- ☒ c: Zero
- ☐ d: one

The current through the solution takes place through the movement of

- ☒ a: Electrically charged particles
- ☐ b: Chemically charged particles
- ☐ c: Magnetically charged particles
- ☐ d: None of these

The SI unit of conductivity is


- ☐ a: Ohms
- ☐ b: Ohms per meter
- ☒ c: Siemens per meter
- ☐ d: None of these

Submit Quiz

8 out of 8


Selected Values:

1. pH value: 9.2

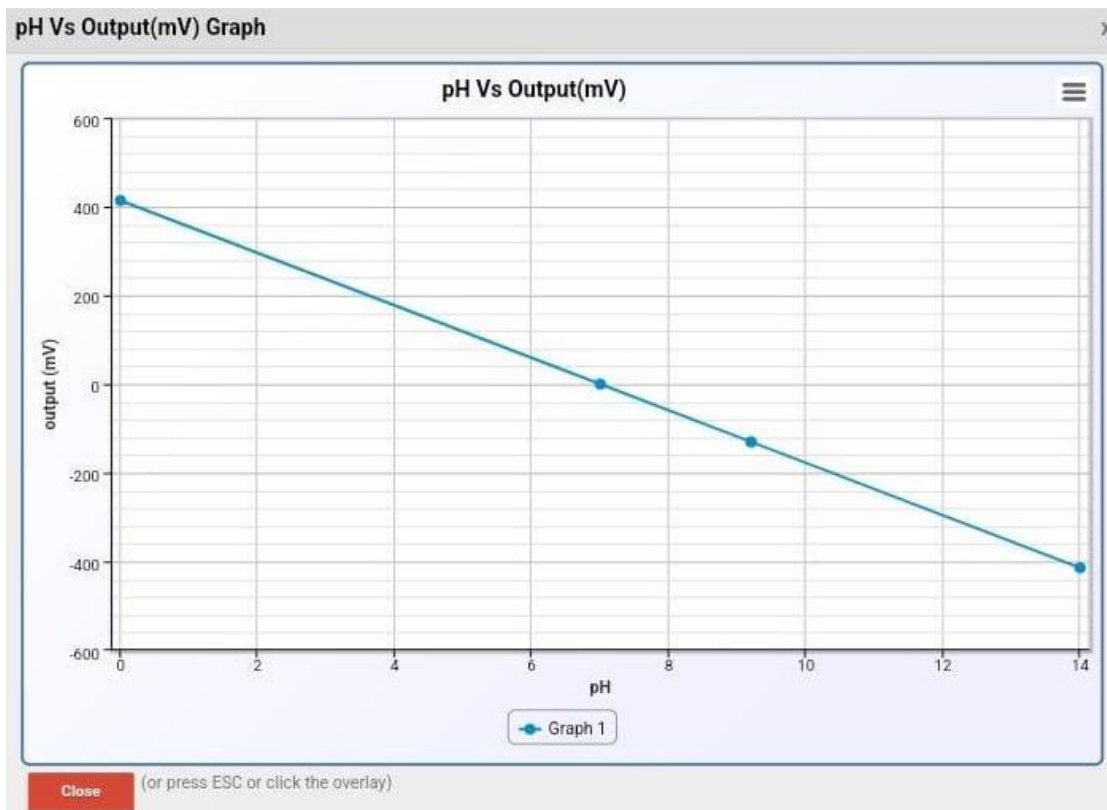


Selected Values:

1. Sample: Sea Water
2. pH Value: 8
3. Output Voltage: -59.16
4. Temperature: 24° C

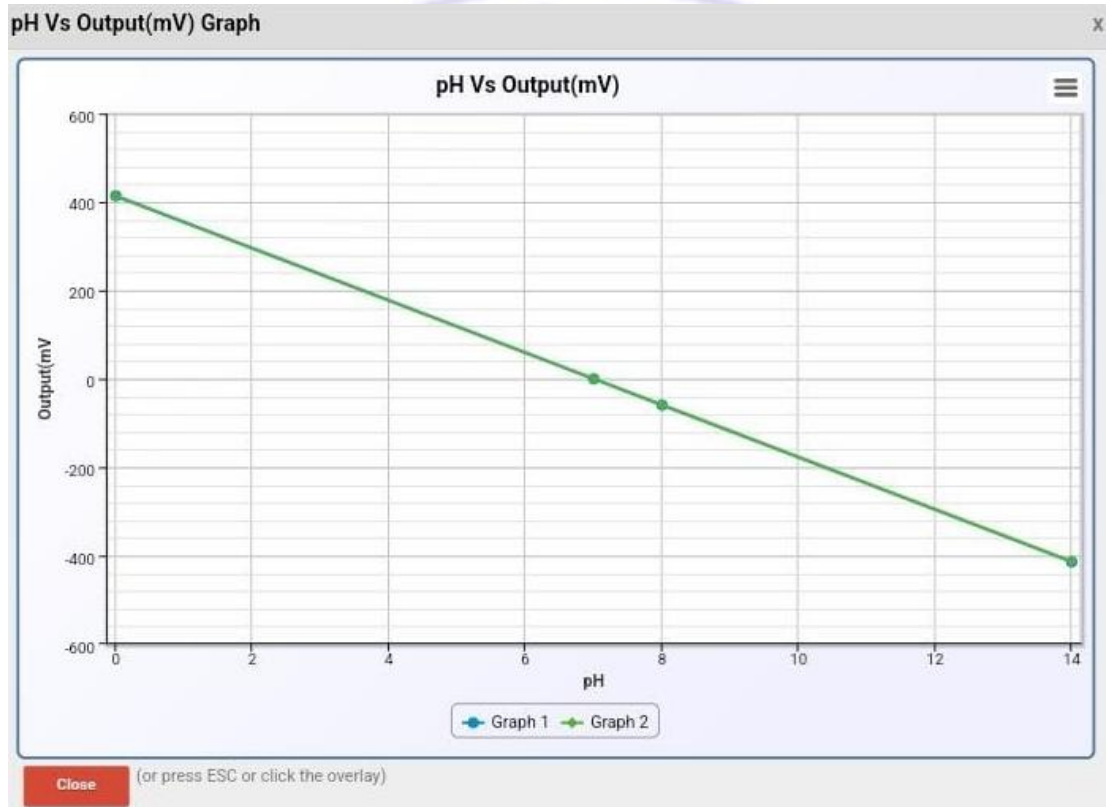


Selected Values:

1. Sample: Orange Juice.
2. pH Value: 3.7
3. Output Voltage: 195.23
4. Temperature: 35° C

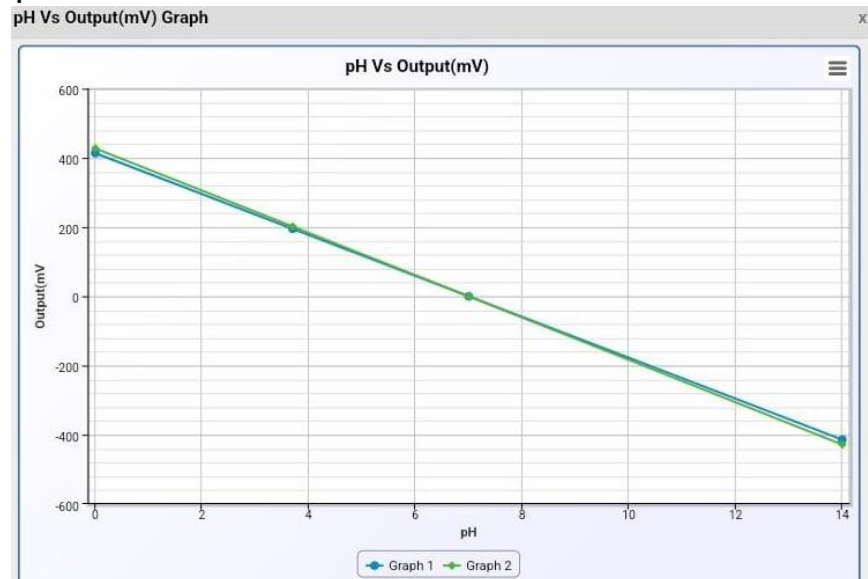

Measurement of Conductivity:

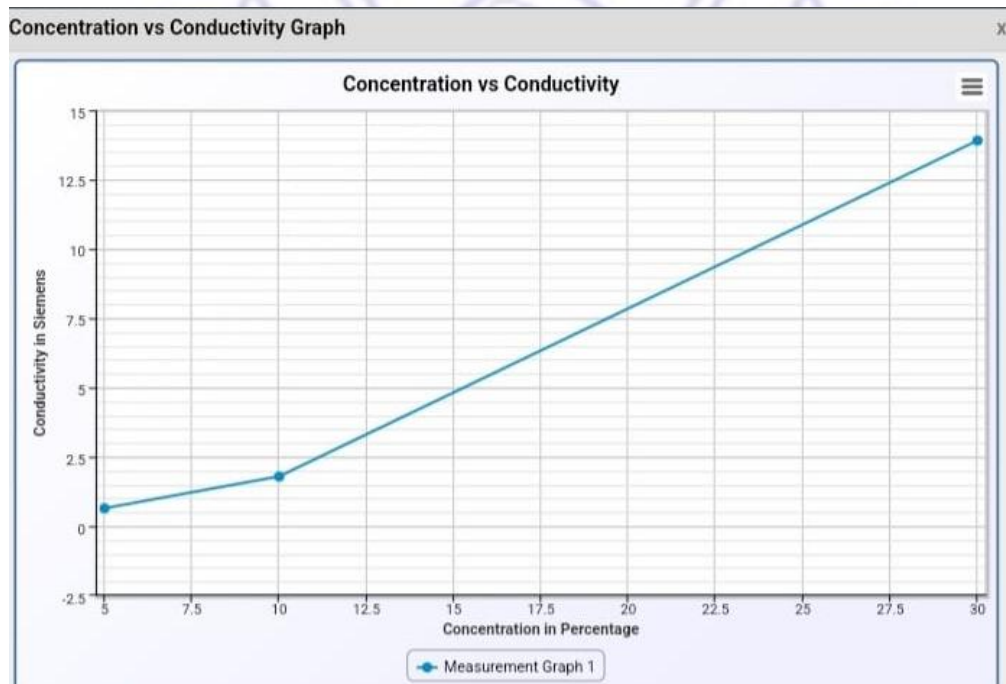


— HCL
Specific conductance at $25^{\circ}\text{C} = 13.91$
Concentration = 30%

1) Temperature = 46°C
18.35 siemens

2) Temp = 20°C
12.85 siemens


3) Temp = 23°C
13.49 siemens

Temperature vs Conductivity Graph

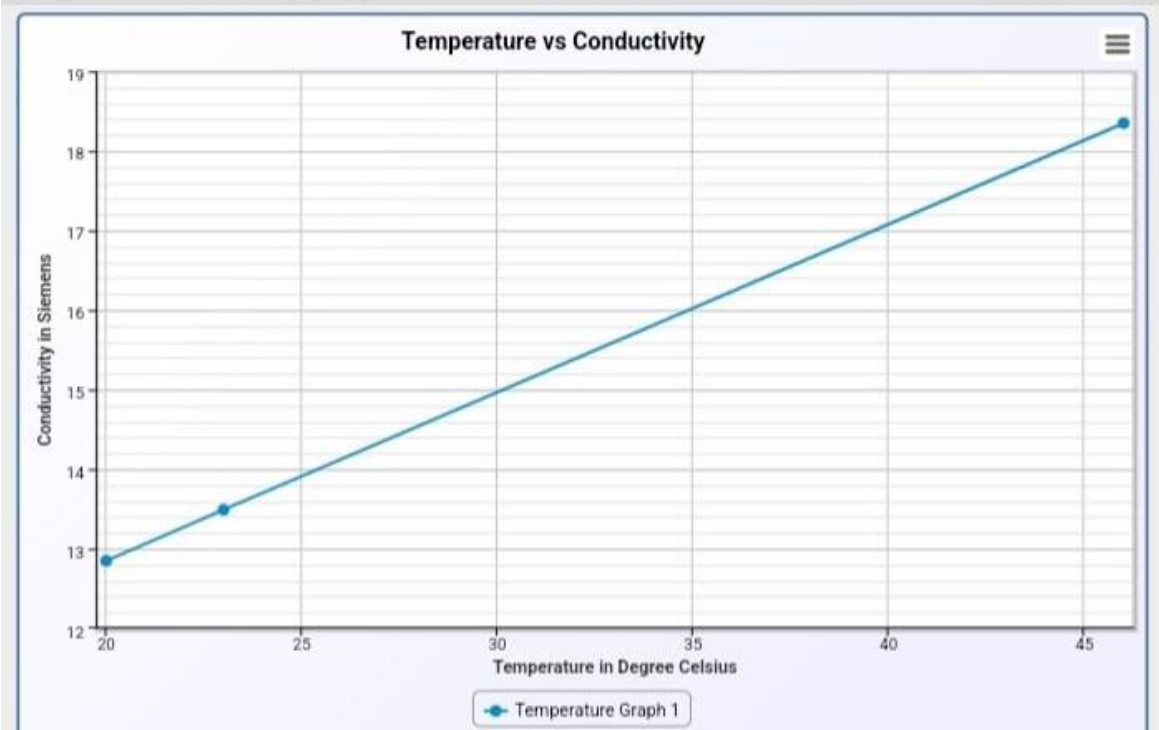

Level-3 Contamination

<--Level-2

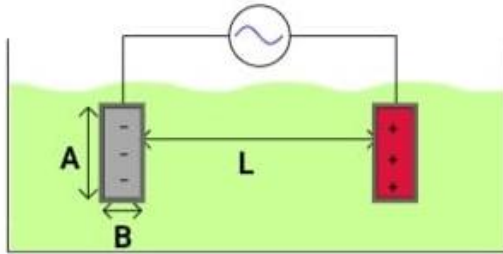

Console window:

Assumption: Half of this value is deposited on each electrode.

Default values for L, A and B:

L=1cm , A=10cm , B=0.1cm

Selected Contamination is 0.3

Control Panel:

Sample:

Concentration:

Cellconstant:

Contamination:

Modified Cell Constant :


Specific Conductance
at 25°C : (Siemens)

Modified specific
conductance value: (Siemens)

[Reload](#)



Posttest:

Virtual Labs
An MoE, Govt of India Initiative

HOMEPARTNERSCONTACT

Electrical Engineering > Sensors Modeling & Simulation Lab > Experiments

Aim

Theory

Pretest

Procedure

Simulation

Posttest

References

Feedback

Simulate the performance of a chemical sensor

pH is a measurement of the activity of -----ions in the sample.
☐ a: Nitrogen ion
☐ b: Carbon ion
☐ c: Oxygen ion
☒ d: Hydrogen ion

In pH measurement Optimum measurement accuracy depends on
☐ a: Solution concentration
☒ b: Salt bridge and Measuring electrode
☐ c: Probe Length
☐ d: All of these

Electrolytes are substances containing
☒ a: Ions
☐ b: cations
☐ c: anions
☐ d: None of these

Electrical conductivity is a measure of the ability of a solution to carry a -----.
☐ a: Inductance
☐ b: Resistance
☒ c: Current
☐ d: Capacitance

In conductivity measurement -----migrate to the negative electrode, the ----- to the positive electrode.
☐ a: Cations,Cations
☒ b: Cations,Anions
☐ c: Anions, Cations
☐ d: None of these

The accuracy of conductivity measurements can be influenced by
☐ a: Polarization
☐ b: Contamination
☒ c: Both, Polarization and Contamination
☐ d: None of these

Submit Quiz

6 out of 6





Community Links

Contact Us

Follow Us

Sakshat Portal
Outreach Portal
FAQ: Virtual Labs

Phone: General Information: 011-26582050
Email: support@vllabs.ac.in



AGPL 3.0 & Creative Commons (CC BY-NC-SA 4.0)

Conclusion:

We studied the characteristics of pH and Conductivity. We also studied effect of temperature on pH measurement and effect of contamination on conductivity measurement.